
# braintronics®

Body&Mind syncSystem

CASADA<sup>TM</sup>  
healthcare

< timeto relax
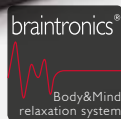
# BETASONIC II

## < timetorelax

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## < timetorelax

Since the company was founded back in the year 2000, the name CASADA has become synonymous with high-quality products in the wellness and fitness sector.

Today, CASADA products are available worldwide. Continuous development, unmistakable design and a strict focus on improving quality of life are just as much characteristic of CASADA products as reliable and top-class quality.

## Congratulations!

By purchasing this massage chair you have demonstrated an awareness for your health.

So that you can benefit from the advantages of this device for a long time, we would ask that you read and observe the safety instructions carefully.

We very much hope you enjoy using your personal BetaSonic II.

## Table of Contents

Safety instructions .....	24
Contraindications.....	25
Equipment   Functions.....	26
Remote control.....	27
Starting up.....	32
Specials.....	35
Troubleshooting .....	39
Storage   Maintenance.....	39
Specifications.....	39
EU – Declaration of Conformity .....	40
Warranty conditions .....	40

# Safety instructions

**Please read the following instructions carefully before using your massage chair, to ensure that it functions properly and with optimal efficiency. Please retain these operating instructions for further use!**

- The massage chair conforms to recognised technical principles and the latest safety regulations.
- The massage chair is equipped with heated sections. Persons who are sensitive to heat should take care when using the product. (ICE 60335-2-32)
- The massage chair does not require maintenance. Any possible repairs may only be carried out by authorised experts.
- Improper use and unauthorised repairs are not permitted for safety reasons and lead to loss of warranty.
- In order to avoid injury please do not stick your fingers between the massage rollers.
- Never touch the mains plug with wet hands.
- Do not allow the massage chair to come into contact with water, high temperatures or direct sunlight.
- Do not use the massage chair in an environment with high humidity levels e.g. in the bathroom.
- Do not use the massage chair in an area with insufficient space or in an area where the ventilation is covered and a supply of fresh air is not guaranteed.
- To avoid the danger of short-circuits, disconnect the massage chair from the power supply if it is not in use for a long period of time.
- Do not use any damaged cables, plugs or loose sockets.
- If cables/plugs are damaged, they must be replaced by the manufacturer, a service representative or by qualified personnel.
- In the case of malfunction, immediately disconnect the product from the mains.
- The company disclaims any liability for damage if the product is in use inappropriately or incorrectly.
- Do not exceed the daily recommended massage duration of 30 minutes in order to avoid excessive strain on the muscles and the nerves.
- To avoid damaging the massage chair, never use pointed or sharp objects.
- In order to reduce the danger of short circuits or fire, it is essential to use only sockets that fit the plug and to ensure that the plug is fully inserted into the socket.
- Disconnect the massage chair from the power supply after use and before cleaning.
- Do not allow children under 14 to clean or maintain the product without supervision.
- If an unexpected power failure occurs during use, immediately set the switch to the OFF position and disconnect the plug in order to avoid damage to the massage chair.
- Do not use the massage chair if the electric cables are wet or damaged.
- Children under the age of 14, people with disabilities and frail people must be supervised at all times when in the vicinity of the massage chair or when using the massage chair.
- Do not use the massage chair immediately after consuming a meal or when under the influence of alcohol.
- Only use the product for the purpose described in the operating instructions.
- Please only use accessories provided by the supplier.
- To prevent damage, please do not stand, sit or bounce on the footrest.
- To avoid damage or injury, do not try to sit on the armrests or backrest.
- Please clean your massage chair regularly and prevent foreign bodies from penetrating the crevices of the massage chair.
- Please do not pull directly on the mains cable.
- When adjusting the position of the footrest or the backrest, please ensure that there are no other objects underneath, in particular no children under the age of 14.
- The product is only intended for domestic use.
- If the flexible power cable becomes damaged, another flexible cable designated for the same purpose must be used or a suitable replacement cable must be purchased from the manufacturer or their service point.

# Contraindications

**If you are unsure whether you should use the massage chair, consult your doctor by all means.**

- Use of the massage chair is prohibited during pregnancy or in the case of one or several of the following complaints: Recent injuries, thrombotic diseases, any kinds of inflammations and swellings, and cancer. If you are undergoing treatment for illnesses or conditions it is recommended that you consult with a doctor before getting a massage.
- Persons with physical, cognitive or mental impairments or persons lacking experience and knowledge should only use this chair under the supervision and guidance of persons responsible for their safety.
- In order to avoid any discomfort or danger when using the massage chair, we would ask the following people to consult a doctor before using the product:
  - Non-ambulant persons
  - Persons receiving medical treatment
  - Patients with back pain and persons with spinal problems / curvature of the spine.
  - Persons with pacemakers
- If you experience pain during the massage, stop using the product immediately.
- If you injure yourself whilst using the chair, seek immediate medical treatment.
- Please do not massage any reddened, swollen or inflamed areas of the skin.
- This massage chair is intended to be used only for home wellness massages and is not a medical massage device. It can therefore not be used as a substitute for specialist medical treatment.
- We recommended that you use the massage chair for 30 minutes on a daily basis. Individual areas of the body should be massaged for a maximum of 15 minutes at a time, as excessive massaging can otherwise be counter-productive.
- This product is not suitable for children under the age of 14.
- Children under the age of 14 must be supervised. Ensure that children do not play with the massage chair.



## “BRAINTRONICS” CONTRAINDICATIONS

- Psychoses (e.g. Schizophrenia, bipolar disorder, endogenous depression)
- Personality disorders
- Epilepsy and similar seizure disorders
- Heart conditions
- Diseases of the central nervous system
- Thrombosis
- Endogenous depression
- ADHD (under-active)
- Recent heart attack or stroke
- Mental disability
- Addictions such as to drugs or alcohol or dependency on medication
- Regular intake of psychotropic drugs
- Pregnancy

**NOTE:** For ethical and legal reasons, braintronics should not be used on children and young people without the explicit permission of the legal representatives. Persons for whom meditation may conflict with religious beliefs should refrain from using the product.

# Equipment | Functions

1. Pillow
2. Shoulder airbag
3. Pocket for remote control
4. Forearm airbag
5. Seat cushion
6. Calves airbag
7. Foot airbag
8. Backrest cushion
9. Armrest
10. Side cover

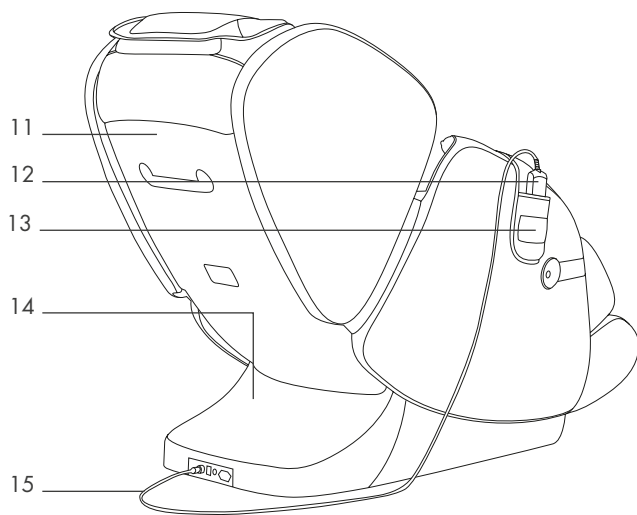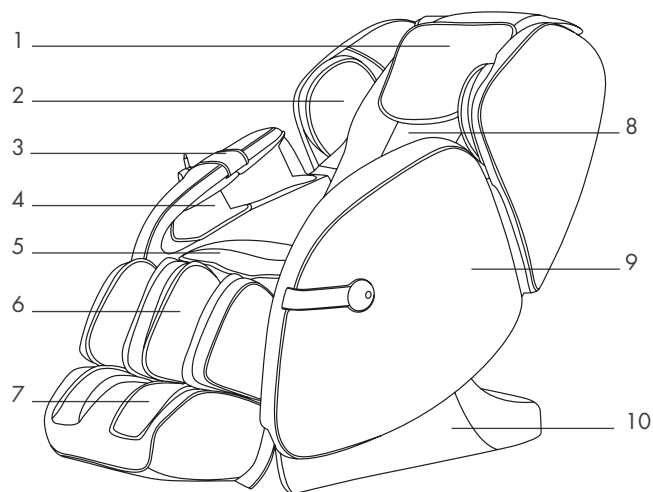
11. Back cover
12. Remote control
13. Mobile phone pocket
14. Safety cover
15. Remote control cable
16. Socket for remote control
17. Power switch
18. Fuse box
19. Power cord connection


# Remote control

1. ON | OFF
2. Menu
3. Auto Programme
4. Manual setting
5. Move the footrest up | down
6. Weightless position (Zero G)
7. Pause
8. Back
9. Advanced setting
10. Move the backrest and footrest up | down


- SCREEN SETTING
- Shoulder height
- Setting



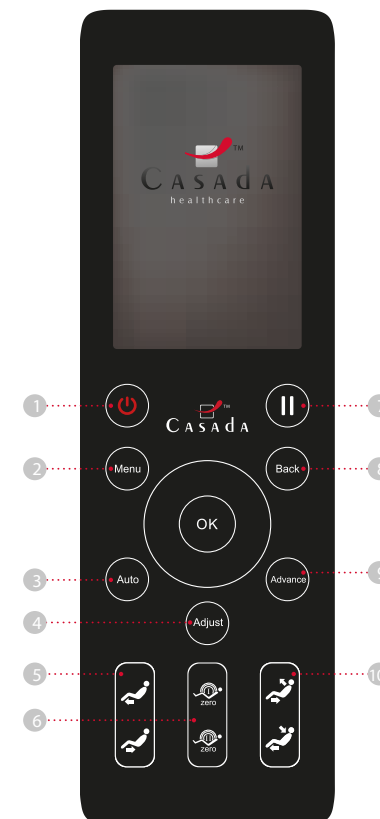
# Remote control


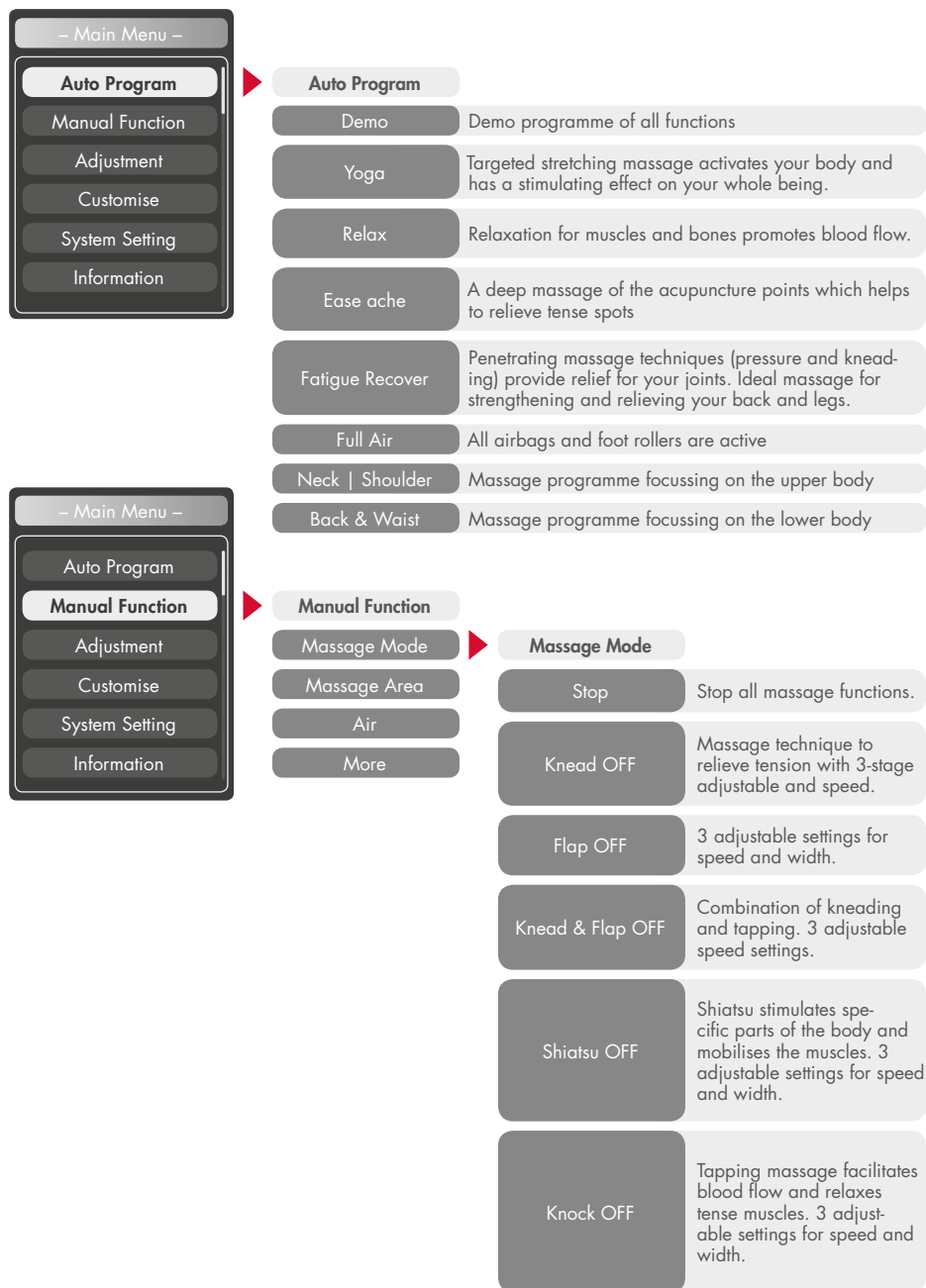
# Remote control
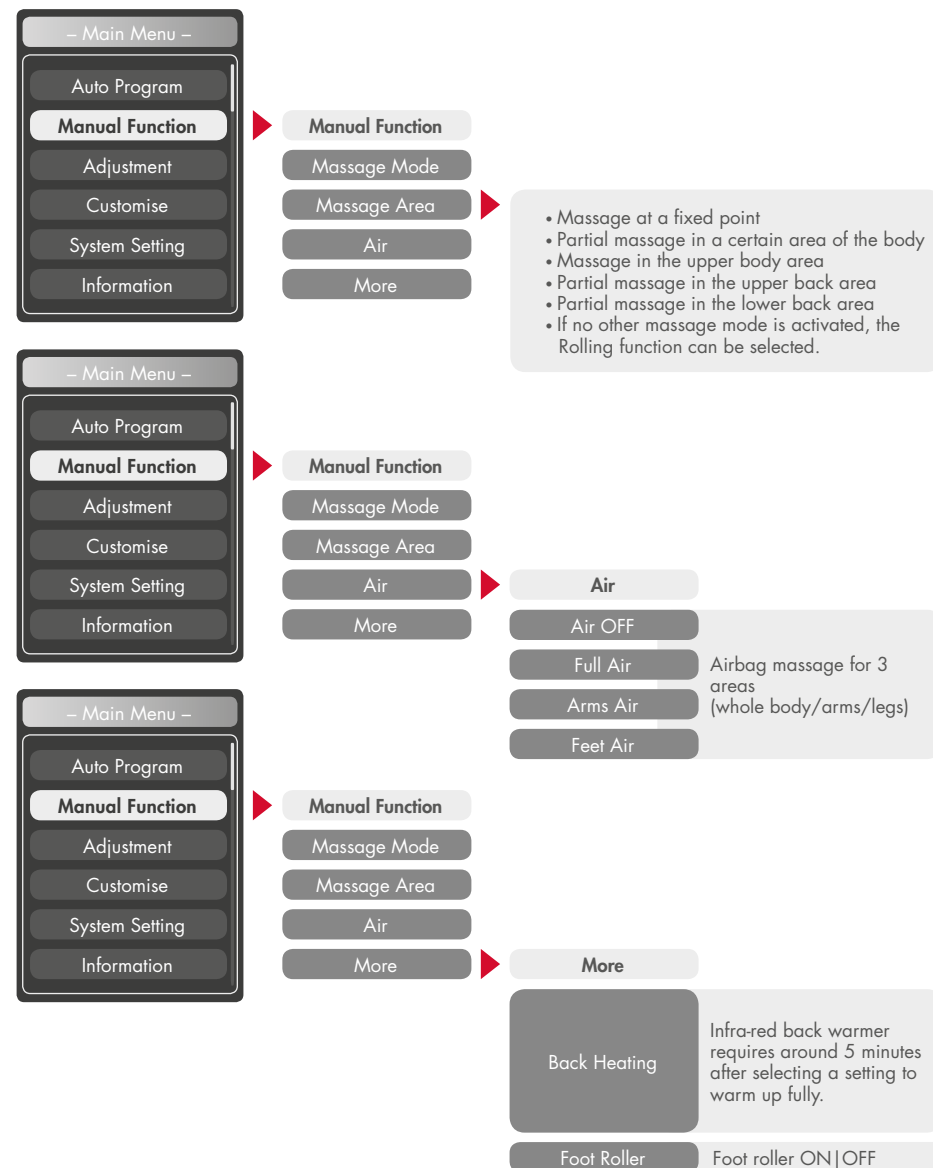
# Remote control
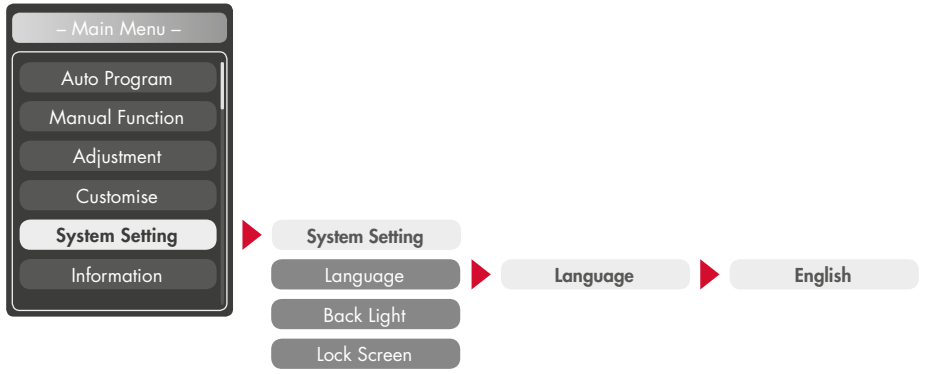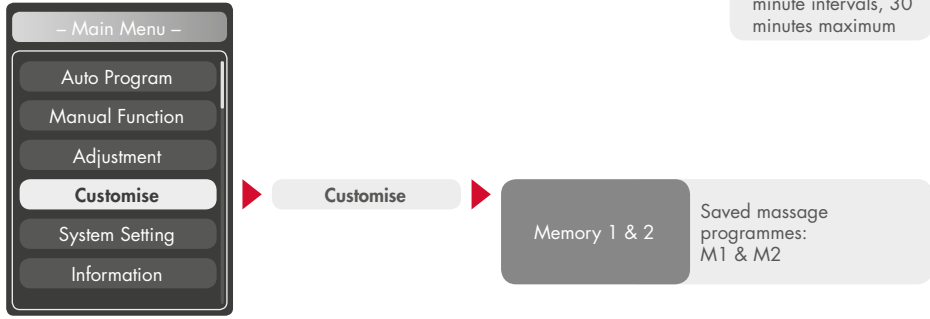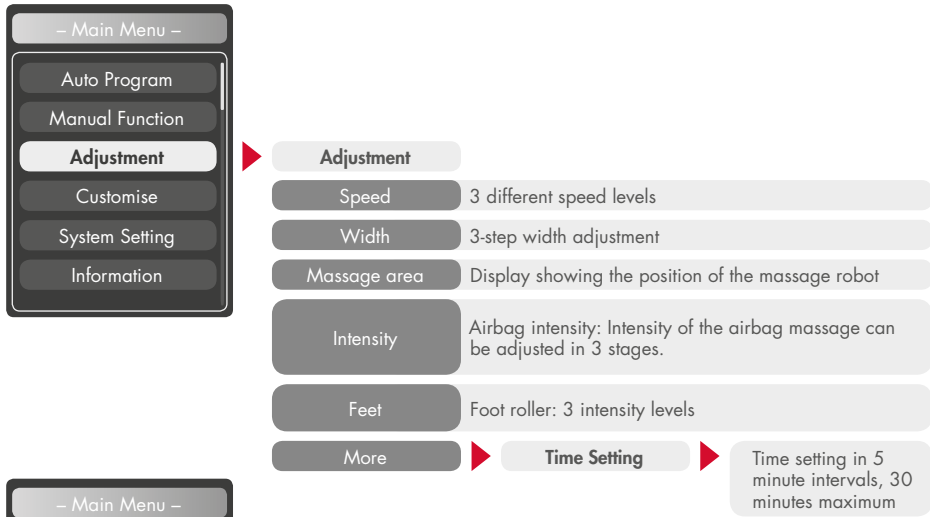
# Remote control
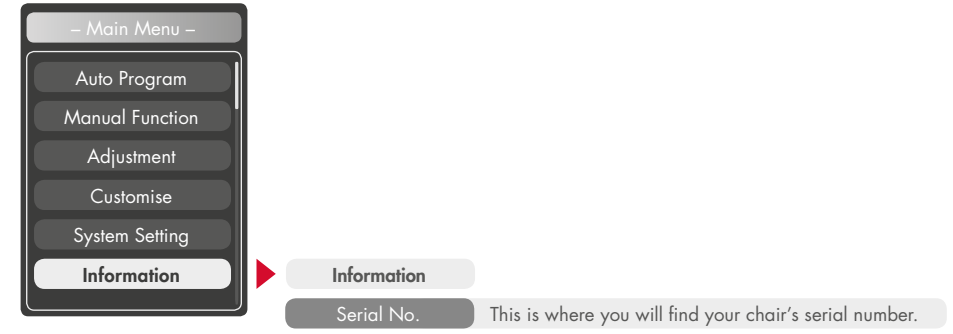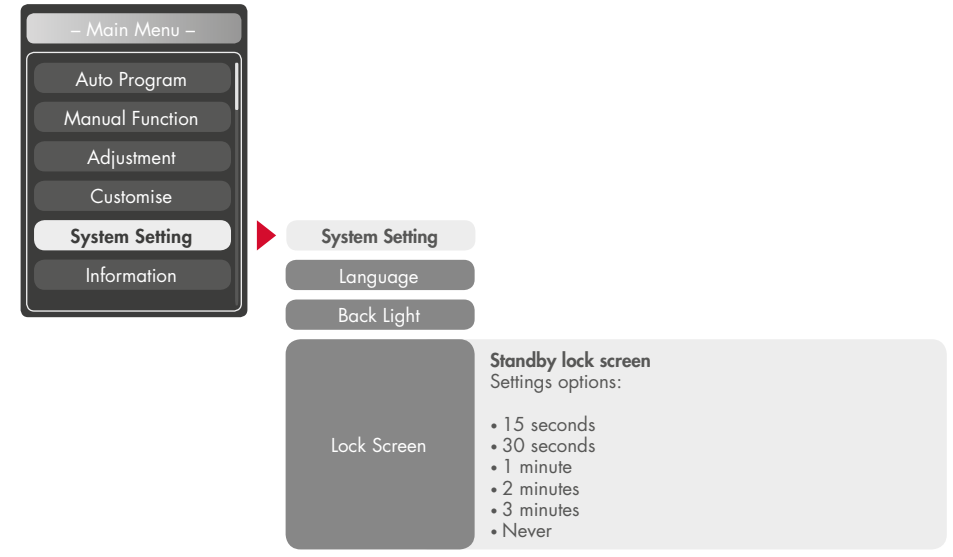
# Starting up

## 1 – ATTACHING THE FOOTREST

- a) First attach the rotating shaft to the fixed axle. Ensure that the R-bolts are the correct way round (Illustration 1).
- b) Once the fixed axle is aligned and combined (Illustration 2), join the pieces together.

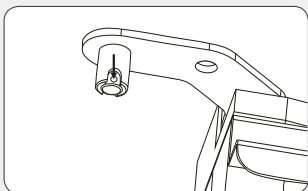
Illustration 1

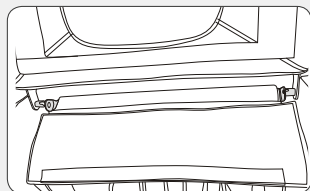
Illustration 2

- c) Place the R-bolt through the holes on both sides to secure the footrest. (Illustration 3)
- d) Assembly complete (Illustration 4)

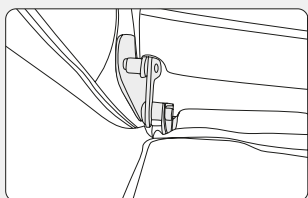
Illustration 3

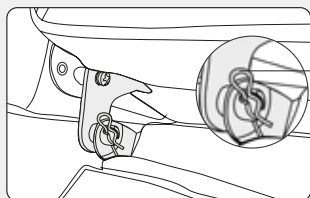
Illustration 4

- e) Connect the air tube and the power connection holder of the foot section to the counterpart coming from the lower seat section.

### NOTE

- To avoid injuries and damage to the product, the footrest should not be lifted from the outside or pressed down. (Illustration 1)
- In order to avoid sitting in an uncomfortable position, make sure that the massage mechanism has not stopped in an upright position in the seat cushion. (Illustration 2)

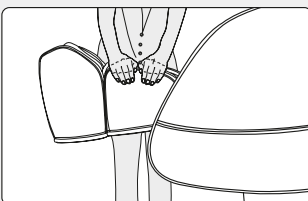
Illustration 1


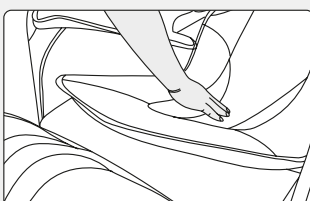
Illustration 2

# Starting up


## 2 – FITTING AND USING THE HEAD AND BACKREST CUSHION

Please connect the backrest cushion as shown in (Illustration 1) using the practical zipper on the massage chair. Please fasten the pillow using the Velcro strips on the upper part of the back cushion (Illustration 2). To clean the backrest cushion, you can easily separate it from the chair by pulling the zipper. We would strongly advise that you do not use the massage chair without the backrest cushion.
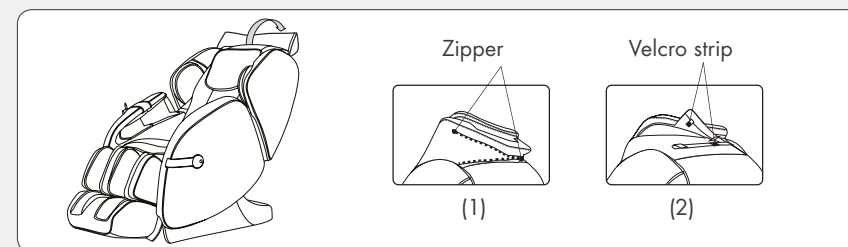
## 3 – FLOOR PROTECTION

In order to prevent excessive wear and tear of your floor covering due to the heavy weight of the massage chair, we recommend placing carpet or similar underneath the chair.
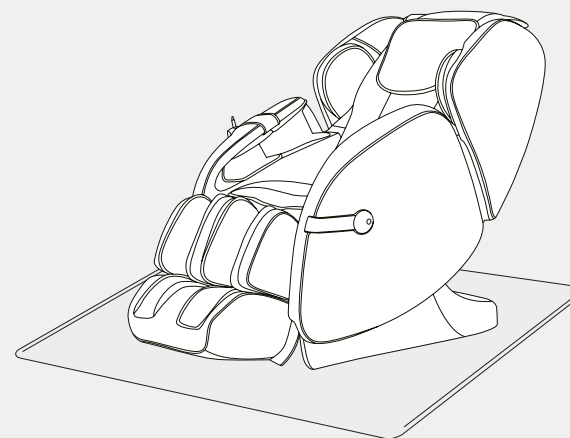
### NOTE

The massage chair should not be pushed over uneven or delicate floor surfaces, but should be carried by at least 2 people. To do this, lift the chair using the handle on the back section and the footrest. This way the footrest is completely folded up.

# Starting up

## 4 – SWITCHING ON THE POWER SUPPLY

Connection procedure

Electrical power switch

## 5 – FINISH MESSAGE

- a) Pressing the power button during a massage stops all massage functions, and both the backrest and the footrest return to their original position.
- b) Unplug the massage chair once you have finished using it.


Switch position and power supply

Connection procedures

## 6 – TRANSPORTATION

Ensure that all cables are far from the ground. Slide the seat on the backrest to the desired position and carefully set it upright again.
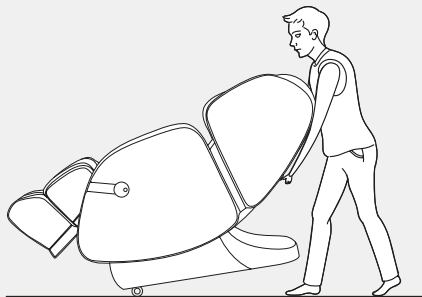
### NOTE

Before transportation the electric power supply must be switched off and the power cable with the remote control cable must be disconnected.

# Specials

## SHIATSU

Shiatsu (finger press) is a form of body therapy originating in Japan. The massage technique consists of soft, rhythmic, far-reaching stretches and rotations. The aim of Shiatsu is to simulate individual parts of the body and to mobilise the muscles.

## 2D MASSAGE

The intelligent 2D massage mechanism allows the massage chair to measure your back to determine a more exact position of the sections of the spine.

## KNEADING

Skin and muscles are either clasped and kneaded between the thumb and the index finger or with both hands. This massage technique is used in particular to ease tension.

## FLAP

Circulation in the muscles is stimulated and the tendons are encouraged to stretch by means of light tapping with a cupped hand, which also reduces tension at the same time.

## TAPPING

Short, tapping movements are performed with the edge of the hand, the palm or the fist. This promotes blood circulation in the skin and softens tense muscles. If tapping is performed at lung height, this can improve the release of mucus in the lung.


## KNEAD & FLAP

A massage combination of kneading and flapping.

## ARTHROSE TRACTION

Arthrosis Traction consists of gentle stretching pressure. It has a stimulating effect on the cartilages, which enhances the transport of bodily fluids, nutrients absorption and can alleviate symptoms of arthrosis.
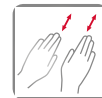
## ROLLING

Pleasant alternation between relaxation and releasing the muscles – also perfectly suited for a final massage, for relaxation and recovery of the muscles.

## CARBON HEATING

Deeply penetrating Carbon infrared heat, blood vessels expand, the blood flow increases and the muscle tone is reduced. In addition, it carries out a highly beneficial effect on nerve cells, which are responsible for reporting pain to the brain.

# Specials

## GUA SHA

The scraping technique of -GuaSha- with gentle pressure opens up superficial blood vessels and stimulates blood circulation. GuaSha is often applied in combination with acupressure 'shu points' for additional stimulation of the whole body for an extra soothing massage experience.

## AIR COMPRESSION

Rising and receding airbags have a pumping motion which produces natural muscle reflexes. This is particularly beneficial in the calf and foot area.

## BLUETOOTH

Wireless data transmission.

## BRAINTRONICS

Braintronics® technology synchronises your brainwaves through audio stimulation. Combined with a specially developed massage programme your stress-ridden body will be able to relax. No matter whether you'd like to relax, dream or learn – anything's possible.

## ADJUSTABLE PILLOW

If necessary, the double cushion can be placed to reduce the neck and shoulder massage intensity.

## SL-SHAPE

Thanks to the „SL“ curve rail technology, which has been approved from human engineering, the range of the massage technology is up to 135cm. This method allows you to enjoy a personal full-body massage from head to thigh and supports the maintenance of the „S“-shaped spine.

## ZERO SPACE

For the lying-position the massage chair moves forward and therefore requires very little distance to the wall behind. Space saving and comfortable!

## ZERO GRAVITY

You put yourself in the position of an astronaut – that "weightless feeling".

## AUTOMATIC

9 specifically developed automatic programmes for the whole body, which offer a big range of different massage techniques. Additionally you can also turn on the heating function, if you'd like.

# Specials

## MEMORY

Using the intelligent memory function the massage chair remembers your favourite massage exactly and begins immediately without a scan. You have the opportunity to save two programmes with different buttons (M1 & M2).

## BACK AREA

An individual applicable to any size massage area.


## ADJUSTABLE FOOTREST

The footrest can be extended, this way it is adaptable to different body sizes.


## RECLINING POSITION

The adjustability of the backrest and footrest contributes to personal well-being of the body position.
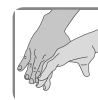
## HUMAN HANDS FEELING

Specially designed massage heads imitate the hand motion of a professional masseur. The maximum stretch space is 6,5 cm, and the maximum stretch angle is 41 degrees, so the mechanical hands can do effective massage to human cervical vertebra thoracic vertebra and lumbar vertebra.
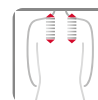
## ADJUST SHOULDER

The massage heads can be moved upward or downward to a suitable position in the shoulder area.

## AIR INTENSITY

Adjustable airbag massage intensity.

## AIR AREA

Three air pressure modes (full body, arm and leg / lower body).

## SPOT


The massage heads can massage a particular spot.

## SPEED

Massage speed is adjustable.

# Specials
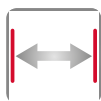
## WIDTH

The distance between the both massage heads can be set (narrow, medium and wide).

## ACUPRESSURE POINTS

Key acupressure points are massaged.

## YOGA

Targeted stretch massages activate the body and have an invigorating effect on the whole of the body.

# Troubleshooting

## FAULT

- The noise of the engine is audible during use.
- The remote control does not work properly.
- The device stops operating.
- The braintronics programmes won't start.

## EXPLANATION

- This is a normal operating noise.
- Ensure that the mains plug is properly inserted into the socket and that the power switch is switched on. Please check the cable of the remote control.
- The chair switches itself off automatically once the pre-set runtime has been reached or after excessive usage. This ensures the chair has a long operational life.
- Make sure that the braintronics app is installed and that your account is activated. In addition, please check whether the Bluetooth® connection between your smartphone and the massagechair persists.

# Storage | Maintenance

- If the device is not used for a longer period of time, the power cable should be wound up and the device should be stored in a dry environment.
- Never store the device at high temperatures or in proximity to open fire. Avoid prolonged exposure to sunlight.
- Please only use commercially-available artificial leather care products for cleaning. Do not clean it with thinner, benzene, or alcohol at any time.
- If the device is on an uneven surface, do not roll or push it; lift it up in order to move it.
- Please take breaks after each use and do not use the device a number of items without pauses.

# Specifications

Dimensions: upright: 76 x 110 x 143 cm  
reclining: 76 x 98 x 168 cm  
Weight: 92 kg  
Voltage: 220-240V ~ 50/60Hz  
Nominal power: 150 W  
Runtime: 20 minutes  
Braintronics: 21 minutes

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